

In-Floor Scale

Model SRV947IFS

The SRV947IFS is a flush mounted platform scale. This is a great alternative when floor space is at a premium. The flush surface eliminates tripping hazards and hallway accessibility issues.

The scale's durable frame is coated with an aggregate infused polyurea coating that provides a slip-free surface for animals, reducing the fear that animals often experience with a scale.

The large capacity allows pets to be easily weighed in kennels by taring out the kennel weight.

A built in surface ring makes cleaning underneath the scale quick and easy.

CAPACITY

400 Lbs or 182 Kgs

LARGE WEIGHING SURFACE

The scale's large 24" x 48" anti-skid surface makes it ideal for gathering weight on large or small animals.

ALWAYS ON

AC wall power means the scale is always on and ready to weigh.

AUTO-HOLD

Scale will automatically lock in on weight above 1 Kilogram.

WIRELESS CAPABILITY *OPTIONAL*

Scale comes with bluetooth wireless output that communicates with the SR Scales app on android devices.

MEMORY RECALL

Recall the last stored weight.



In-Floor Scale

Model SRV947IFS

SR Scales[®]
by SR Instruments, Inc.


Quality. Reliability. Accuracy.

SYSTEM SPECIFICATIONS


MAX WEIGHT	400 Lb / 182 Kg
DISPLAY RESOLUTION	.2 Lb / .1 kg
BUTTON FUNCTIONS	Zero/Weigh, lb/kg, Hold/Recall
DISPLAY UNITS	Pounds / Kilograms
POWER	120 VAC 50-60 Hz
HELD WEIGHT	Press "Hold/Recall" to store the current weight. Pressing "Hold/Recall" again will retrieve the last stored weight.
ALWAYS ON	Scale is always on in the background zero-tracking and will automatically wake up and lock a stable weight when weight is detected.
WARRANTY	Two (2) Years

ACCESSORIES AND OPTIONS


DESCRIPTION	PART NUMBER
IN FLOOR SCALE WITH WIRELESS	SRV947IFS-BT

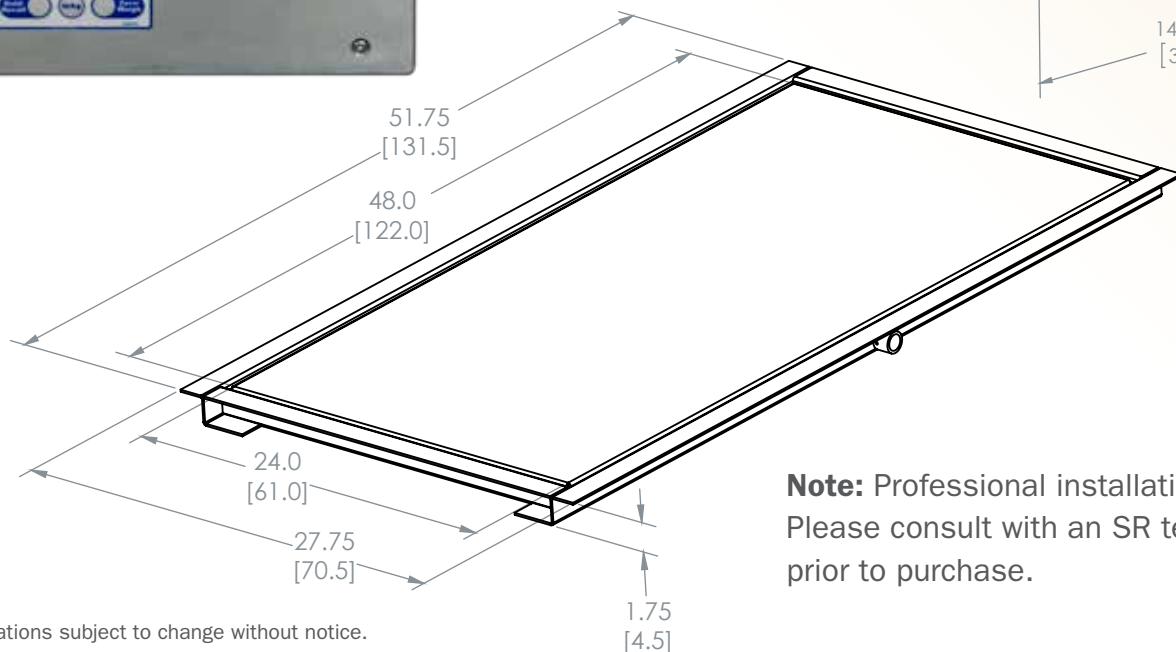

PACKAGING

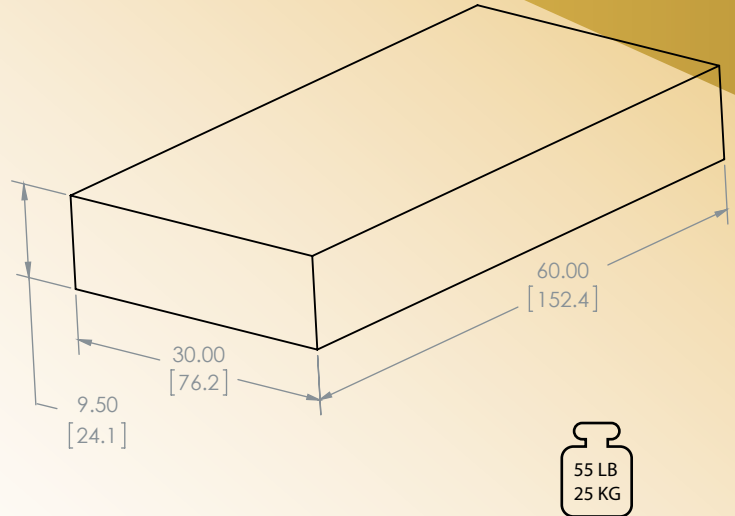

DIMENSIONS in [cm]

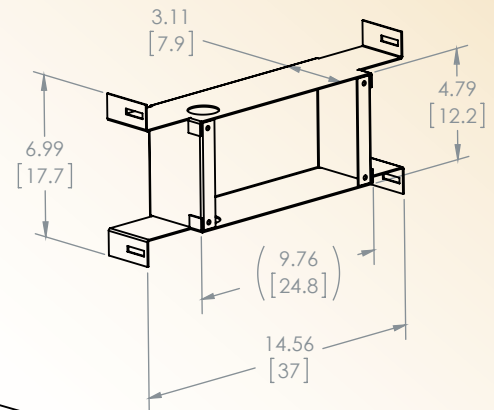

Note: Professional installation is required. Please consult with an SR team member prior to purchase.

Specifications subject to change without notice.

SR Instruments, Inc.

600 Young Street, Tonawanda, New York 14150 USA
T: 716.693.5977 • F: 716.693.5854 • TF: 1.800.654.6360

www.srscales.com

email: info-request@srinstruments.com